



I.D. Technologies

Division of The Hitt Companies, Inc.

MANUFACTURERS OF FINGERPRINTING & IDENTIFICATION PRODUCTS



The Forensic Experts

**FORENSIC
EVIDENCE COLLECTION
FINGERPRINTING SUPPLIES**

**Eliminate the Middleman
Buy Direct From The Manufacturer**

www.idtechnologies.com

Specimen Collection

- 5 Sexual Assault Evidence Collection Kits**
- 9 Drug Facilitated Toxicology Kits**
- 11 Child Physical Abuse Kit**

- 12 DNA Specimen Collection Kits**
- 17 Child Identification**
- 20 DNA FTA[®] Matrix Cards**

- 25 Blood and Urine Specimen Collection Kits**
- 29 Blood Specimen Collection Tubes**
- 32 On-Site Drug Screening Kits**

- 35 Postmortem Kits**
- 36 Postmortem DNA Blood Collection Kit**
- 39 Postmortem Fingerprinting Kits**

- 41 Gunshot Residue Evidence Collection Kits**
- 44 GSR Stub Information**
- 45 Gunshot Residue Collection Supplies**

OUR BEST-SELLING VICTIM SEXUAL ASSAULT EVIDENCE COLLECTION KITS

AVAILABLE WITH EITHER BLOOD TUBES OR A FILTER PAPER BLOOD COLLECTION CARD

Used by agencies throughout the world, I.D. TECHNOLOGIES' Victim Sexual Assault Evidence Collection Kits were developed to provide consistent, thorough evidence in cases of suspected sexual assault. They are available with either blood tubes in the RE-1MS or with a filter-paper blood collection card in the RE-1MS(FS), which eliminates the need to keep the kit refrigerated after specimen collection.



FS-RE1MS & FS-RE1MS(FS)

KIT COMPONENTS

- Kit Box (10¼" L x 6¾" W x 2⅜" H) pre-sealed at I.D. TECHNOLOGIES via box-top integrity sealing tabs
- Kit Instruction Sheet printed front and back

Step 1	Authorization for Collection and Release of Evidence and Information Form printed on 3-ply NCR paper
Step 2	Victim's Medical History and Assault Information Form printed on 3-ply NCR paper
Step 3	Foreign Material Bag containing one 20" x 30" white paper sheet
Step 3	Underpants Bag
Step 3	Outer Clothing Bags (2 ea.)
Step 4	Debris Collection Envelope containing one paper bundle, one package of sterile, cotton-tipped swabs (2/pkg), and one fingernail scraper
Step 5	Pubic Hair Combing Envelope containing one swabbie towel and one 5" black plastic comb
Step 6	Pulled Pubic Hairs Envelope
Step 7	Vaginal Swabs and Smears Envelope containing one plastic slide mailer with two microscope slides, two swab boxes, and two packages of sterile, cotton-tipped swabs (2/pkg)
Step 8	Rectal Swabs and Smear Envelope containing one plastic slide mailer with one microscope slide, one swab box, and one package of sterile, cotton-tipped swabs (2/pkg)
Step 9	Oral Swabs and Smear Envelope containing one plastic slide mailer with one microscope slide, one swab box, and one package of sterile, cotton-tipped swabs (2/pkg)
Step 10	Pulled Head Hairs Envelope
Step 11	Known Saliva Sample Envelope containing one stamped and folded filter paper disk
Step 12	FS-RE1MS - Known Blood Samples Envelope containing one 7 ml ACD (yellow-stoppered) and one 7 ml EDTA (lavender-stoppered) blood collection tube in a bubble pack blood tube protector
	FS-RE1MS(FS) - Known Blood Samples Envelope containing one filter paper blood collection card, one alcohol prep pad, one sterile gauze pad, one 21 gauge fixed-depth sterile lancet, and one adhesive bandage
Step 13	Anatomical Drawings printed on 3-ply NCR paper
Also included: Police Evidence Seals for resealing kit box after evidence collection, a Law Enforcement Forms Envelope on base of kit box, one Biohazard Label, and an FDA Product Insert	

Item.	Description
FS-RE1MS	Victim Sexual Assault Evidence Collection Kit (blood collection tubes), 13 kits/case
FS-RE1MS(FS)	Victim Sexual Assault Evidence Collection Kit (blood collection card), 13 kits/case

FS-RE1MS - Victim Sexual Assault Evidence Collection Kit - Blood Tubes



FS-RE1MS(FS) - Victim Sexual Assault Evidence Collection Kit - Finger Stick



Tip

In cases where a drug facilitated sexual assault is suspected, use of our companion DFRE Kits on the following page may be appropriate. Designed to be used in conjunction with the RE-1MS, RE-1MS(FS), and many other sexual assault kits, the DFRE Kits will help you collect evidence in a proper manner and provide appropriate chain of custody procedures.

DRUG FACILITATED SEXUAL ASSAULT TOXICOLOGY KITS

The use of "Date-Rape Drugs" such as Rohypnol™ ("Ruffies") and Gamma-Hydroxybutyrate ("GHB") is being reported more frequently in sexual assault cases throughout the country. Because of the amnesia effect of these drugs, the victim may not be able to remember all the details of the assault or may have difficulty identifying the assailant(s).

Our Drug Facilitated Sexual Assault Specimen Collection Kits (DFRE) are designed to be used in conjunction with your existing sexual assault evidence collection kits when victimizing drugs are suspected. The kit instruction sheets are written for medical personnel and assure proper collection and chain of custody. After the examination the kit is handed over to the investigating officer, along with the sexual assault evidence collection kit, for submission and analysis by the crime laboratory.

Blood and Urine Specimen Collection Kit

FS-BUDFRE	
KIT COMPONENTS	
1 ea.	Pre-sealed Kit Box (6" L x 4¾" W x 2½" H) with chain of custody and biohazard label
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Victim Information Form, printed on 2-ply NCR paper
1 ea.	Inner Box Blood Tube and Urine Bottle Holder
1 ea.	100 ml Urine Specimen Bottle
1 ea.	Multi-sample Needle (20 g x 1½")
1 ea.	Safety Tube and Needle Holder
1 ea.	PVP Prep Pad (alcohol-free)
2 ea.	10 ml Blood Collection Tubes (gray-stoppered)
3 ea.	Tamper-evident Specimen Security Seals
1 ea.	Ziplock Bag containing a liquid absorbing sheet
1 ea.	Tamper-evident Kit Box Shipping Seal



Item.	Description
FS-BUDFRE	Drug Facilitated Blood and Urine Specimen Collection Kit, 10 kits/case

Urine Specimen Collection Kit

FS-UCDFRE	
KIT COMPONENTS	
1 ea.	Pre-sealed Kit Box (5" L x 3¾" W x 2½" H) with chain of custody and biohazard label
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Victim Information Form, printed on 2-ply NCR paper
1 ea.	Specimen Bottle Holder (polyurethane foam)
1 ea.	100 ml Urine Specimen Bottle
1 ea.	Specimen Bag containing a liquid-absorbing sheet
1 ea.	Tamper-evident Specimen Security Seal
1 ea.	Tamper-evident Kit Box Shipping Seal



Item.	Description
FS-UCDFRE	Drug Facilitated Urine Specimen Collection Kit, 10 kits/case

Hair Specimen Collection Kit

Detection of drugs in drug facilitated sexual assault is limited due to their rapid elimination from the victim's urine; however, hair traps and retains date rape drugs and/or their characteristic metabolites, thereby allowing the detection of a variety of drugs used to facilitate a sexual assault for several months after the assault.


Our Hair Specimen Collection Kit has been designed to assure proper hair collection and ensure proper chain of custody and documentation of the evidence. Our step-by-step instruction booklet makes the collection of hair simple and straightforward, thus providing the best specimen for laboratory testing.



FS-HCDFRE

KIT COMPONENTS

1 ea.	Tamper-evident, pre-sealed Kit Envelope, 6" x 9"
1 ea.	Kit Instruction Booklet
1 ea.	Authorization for Collection and Release of Evidence Form, printed on 3-ply NCR paper
1 ea.	Method A Hair Collection Envelope
1 ea.	Method B Hair Collection Envelope
2 ea.	Folded Paper Bindles
1 ea.	Paper-covered Wire Twist Tie (not shown)
2 ea.	Tamper-evident Evidence Security Seals
1 ea.	Kit Mailing Envelope, 7½" x 10½"
1 ea.	Tamper-evident Kit Mailing Envelope Shipping Seal

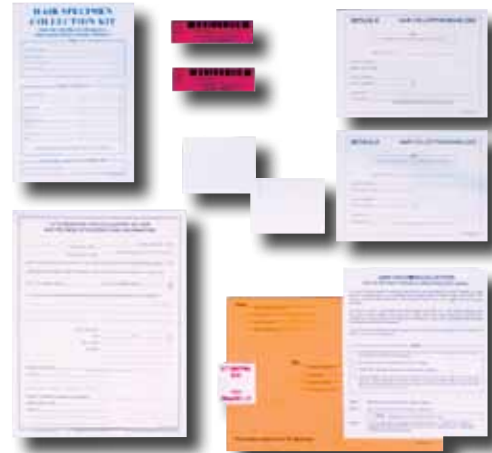

Item.	Description
FS-HCDFRE	Hair Specimen Collection Kit, 10 kits/case

Disposable Swab Drying Rack

Our Disposable Swab Drying Rack is specifically designed for drying the swabs used in Sexual Assault Evidence Collection Kits.

FS-SADR01

FEATURES

	Prevents the possibility of losing track of which swabs are from which body orifice during the air-drying process or during an examination
	Prevents cross-contamination caused by one set of swabs coming in contact with another set during the air-drying process or during an examination
	Holds up to 12 swabs (4 vaginal, 2 anal/rectal, 2 oral, 2 dried secretions, and 2 "other swabs")



Item.	Description
FS-SADR01	Disposable Swab Drying Rack, 50 units/case

SUSPECT SEXUAL ASSAULT EVIDENCE COLLECTION KITS

AVAILABLE WITH EITHER BLOOD TUBES OR A FILTER PAPER BLOOD COLLECTION CARD

Our Suspect Sexual Assault Evidence Collection Kits are supplied with all the necessary components for the collection of specimens from a suspected assailant and are designed to eliminate any confusion as to what specimens should be collected.

The kits are supplied with step-by-step instruction sheets describing the proper collection procedures for suspect apprehension within 24 hours or after 24 hours of the assault.

RS-0GEN - Suspect Sexual Assault Evidence Collection Kit - Blood Tubes



FS-RS0GEN & FS-RS0GEN(FS)

KIT COMPONENTS

- Pre-sealed Kit Envelope, 7½" x 10½"
- Kit Instruction Sheet, printed front and back

Step 1	Two Outer Clothing Bags and One Underpants Bag
Step 2	Foreign Material Collection Envelope containing one paper bindle, one package of sterile, cotton-tipped swabs (2/pkg), and one plastic fingernail scraper
Step 3	Pubic Hair Combing Envelope containing one swabbie towel and one 5" black plastic comb
Step 4	Pulled Pubic Hairs Envelope
Step 5	Pulled Head Hairs Envelope
Step 6	Known Saliva Sample Envelope containing one stamped and folded filter paper disk
Step 7	FS-RS0GEN - Known Blood Samples Ziplock Bag containing one 7 ml ACD (yellow-stoppered) and one 7 ml EDTA (lavender-stoppered) blood collection tube in a bubble pack blood tube protector FS-RS0GEN(FS) - Known Blood Samples Envelope containing one filter paper blood collection card, one alcohol prep pad, one sterile gauze pad, one 21 gauge fixed-depth sterile lancet, and one adhesive bandage
Step 8	Anatomical Drawings Form printed on 8½" x 11" paper
Also included: Police Evidence Seal for resealing kit envelope after evidence collection, a Biohazard Label and an FDA Product Insert	

Item.	Description
FS-RS0GEN	Suspect Sexual Assault Evidence Collection Kit (blood collection tubes), 10 kits/case
FS-RS0GEN(FS)	Suspect Sexual Assault Evidence Collection Kit (blood collection card), 10 kits/case

RS-0GEN(FS) - Suspect Sexual Assault Evidence Collection Kit - Finger Stick




Child Physical Abuse Information and Specimen Collection Kit

This kit was developed to assist medical personnel in the documentation of cases of suspected child physical abuse. The kit contains medical history and physical examination forms as well as anatomical drawings for showing the locations of both new and old trauma on the child's body.

The kit also includes an envelope for photographs, a paper bag for the child's clothing if blood is noted, and a bite/suck mark swab envelope for DNA analysis by the crime laboratory.

FS-CA0GEN	
KIT COMPONENTS	
<ul style="list-style-type: none"> • Pre-sealed Kit Envelope, 9" x 12" • Kit Instruction Sheet with helpful information to assist medical personnel on back 	
Step 1	Authorization for Collection and Release of Evidence and Information Form
Step 2	Four Kit Forms, printed on 2-ply NCR paper
Step 3	Bite/Suck Mark Swab Envelope containing one package of sterile, cotton-tipped swabs (2/pkg)
Step 4	Clothing Bag
Step 5	Photographs/Film Envelope
Also included: Police Evidence Seal for resealing kit envelope after evidence collection and FDA Product Insert	

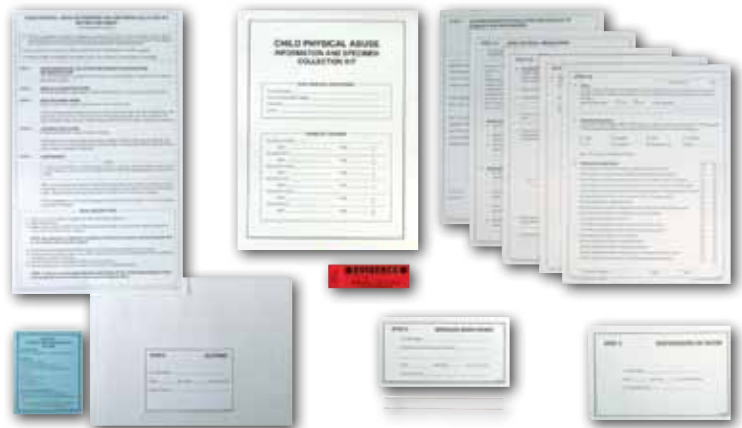

Item.	Description
FS-CA0GEN	Child Physical Abuse Information and Specimen Collection Kit, 10 kits/case

DNA SPECIMEN COLLECTION KITS


With the acceptance of DNA by the courts since the late 1980's, law enforcement agencies rely more and more on DNA evidence in their investigations.

The kits described in this section of our catalog have been specifically designed to assist the investigating officer in the collection and preservation of both "known" and "unknown" DNA specimens.

Based on years of experience designing DNA Database Kits for crime laboratories throughout the country, I.D. TECHNOLOGIES has developed a series of DNA Suspect Specimen Collection Kits specifically for use by investigating officers.

By comparing a suspect's "known" DNA to "unknown" DNA evidence collected at a crime scene, investigating officers can focus their investigation on one suspect rather than several.

Note: Before ordering, check with your crime lab to see what type of suspect DNA specimen best fits their requirements.

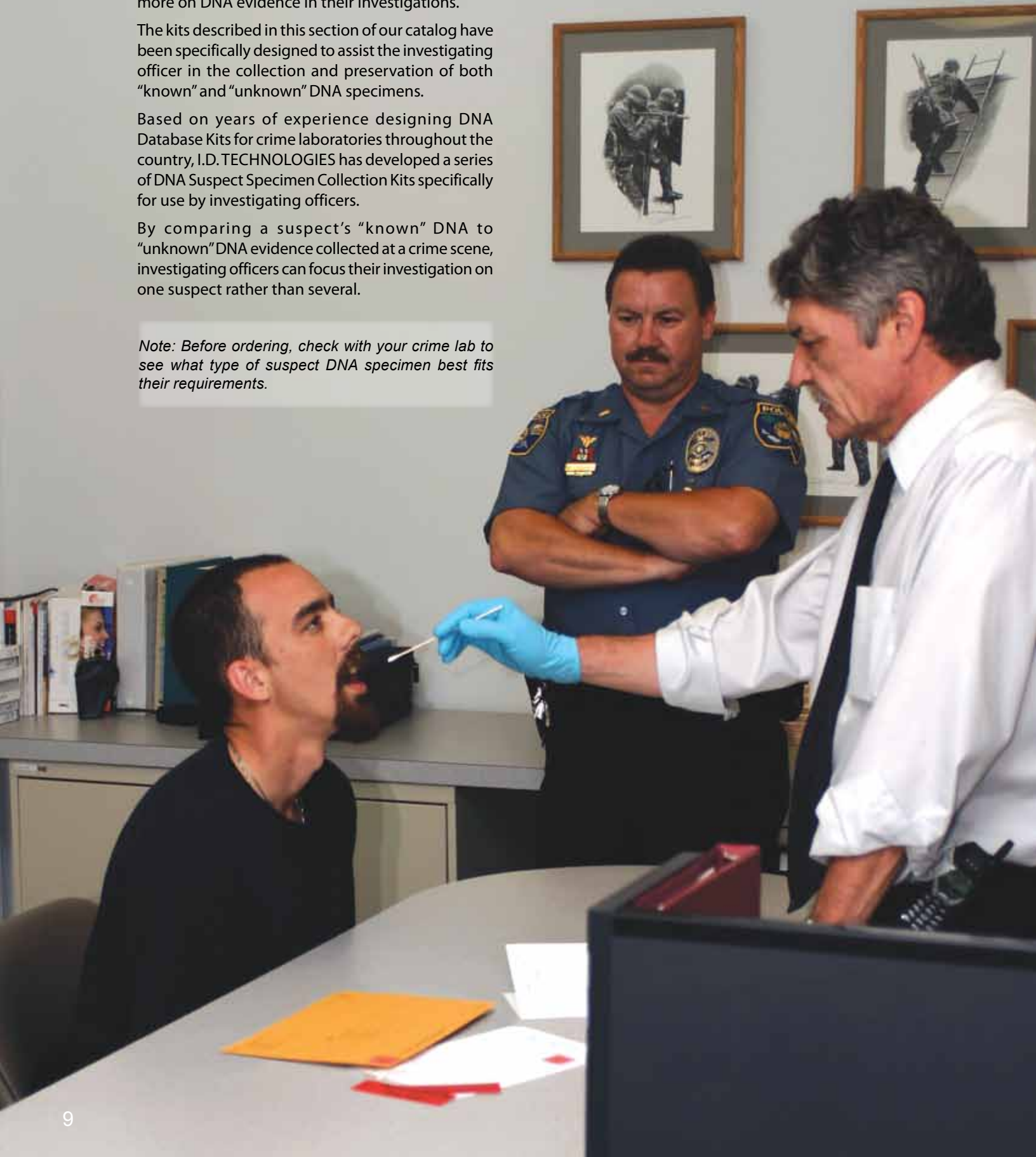

Suspect DNA Buccal Swab Collection Kit



Our Buccal Swab Collection Kit is the least intrusive of our series of Suspect DNA Specimen Collection Kits and is designed for use by nonmedical personnel. The cotton-tipped swabs are placed in the subject's mouth and the inside of the cheek is vigorously swabbed, resulting in the transfer of buccal cells onto the cotton bulb of the swab. The kit is either mailed or hand delivered to the crime laboratory for "known" DNA analysis and comparison with "unknown" DNA from a crime scene.

FS-DNA0GEN(BS)

KIT COMPONENTS

1 ea.	Tamper-evident, pre-sealed Kit Envelope, 6" x 9"
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Buccal Swab Information Card containing a swab envelope
1 ea.	Ziplock Bag containing a desiccant packet
1 ea.	Pair of Nitrile Barrier Gloves
1 ea.	Sterile, cotton-tipped Swabs (2/pkg)
2 ea.	Tamper-evident Police Evidence Seals
1 ea.	Kit Mailing Envelope, 7½" x 10½"
1 ea.	Kit Shipping Seal



Item.	Description
FS-DNA0GEN(BS)	Suspect DNA Buccal Swab Collection Kit, 25 kits/case

Suspect DNA Finger Stick Blood Collection Kit



In most states our Finger Stick Blood Collection Kit can be used by nonmedical personnel. Using the sterile lancet provided, simply stick the subject's left or right ring finger, then squeeze the finger allowing the blood to fill the inner circles on the blood collection card. The kit is then sealed and forwarded to the crime laboratory for analysis.

FS-DNA0GEN(FS)

KIT COMPONENTS

1 ea.	Tamper-evident, pre-sealed Kit Envelope, 6" x 9"
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Fingerstick Blood Information Card containing a blood collection card
1 ea.	Pair of Nitrile Barrier Gloves
1 ea.	Alcohol Prep Pad
1 ea.	Sterile Gauze Pad
1 ea.	21 Gauge Fingerstick Lancet
1 ea.	Adhesive Bandage
1 ea.	Ziplock Bag containing a desiccant packet
1 ea.	Tamper-evident Police Evidence Seal
1 ea.	Kit Mailing Envelope, 7½" x 10½"
1 ea.	Kit Shipping Seal



Item.	Description
FS-DNA0GEN(FS)	Suspect DNA Finger Stick Blood Collection Kit, 25 kits/case

Suspect DNA Whole Blood Collection Kit



Our Whole Blood Specimen Collection Kit requires the blood be drawn by trained medical personnel, however, the responsibilities of the investigating officer and that of the medical personnel are clearly defined and separated to eliminate problems at the time of collection and possibly later in court. After specimen collection, the kit box is sealed and forwarded to the crime laboratory.

FS-DNA0GEN(WB)

KIT COMPONENTS

1 ea.	Tamper-evident, pre-sealed Kit Box (5" L x 3¾" W x 2½" H) with biohazard label affixed
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Blood Collection Information Card
1 ea.	Ziplock Bag containing a liquid-absorbing sheet
1 ea.	7 ml EDTA Blood Collection Tube
1 ea.	Bubble pack Blood Tube Protector
1 ea.	Tamper-evident Blood Tube Security/Identification Seal
1 ea.	Safety Tube and Needle Holder
1 ea.	Multi-sample Needle (20 g x 1½")
1 ea.	Alcohol Prep Pad
1 ea.	Tamper-evident Kit Box Shipping Seal



Item.	Description
FS-DNA0GEN(WB)	Suspect DNA Whole Blood Collection Kit, 25 kits/case

Victim DNA Buccal Swab Collection Kit



Our Victim DNA Buccal Swab Collection Kit is designed for use by nonmedical personnel for the collection of "known" DNA from victims (both adult and minors). Simply swab the inside of the victim's cheek and gum line which will result in the transfer of buccal cells onto the cotton bulbs of the swabs.

After collection the kit is either mailed or hand delivered to the crime laboratory for comparison with "unknown" DNA found at the crime scene to corroborate the victim's testimony.

FS-DNA0GEN(BS)

KIT COMPONENTS

1 ea.	Tamper-evident, pre-sealed Kit Envelope, 6" x 9"
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Victim Buccal Swab Information Card containing a swab envelope
1 ea.	Ziplock Bag containing a desiccant packet
1 ea.	Pair of Nitrile Barrier Gloves
1 ea.	Sterile, cotton-tipped Swabs (2/pkg)
2 ea.	Tamper-evident Police Evidence Seals
1 ea.	Kit Mailing Envelope, 7½" x 10½"



Item.	Description
FS-DNA0GEN(BS)	Victim DNA Buccal Swab Collection Kit, 25 kits/case

Criminal Paternity DNA Specimen Collection Kit



For the past several years, there has been a growing concern with young adults and pedophiles engaging in consensual sexual activities with minors. These relations often lead to the impregnation of the victim. These crimes of statutory rape or incest are often committed by persons in a position of authority or trust. The assailants use their position of power to cover up their crimes by coercing, intimidating, or manipulating their underage victims into undergoing a legal abortion, therefore eliminating any evidence that, through DNA typing, would link them to the minor.

Our Criminal Paternity DNA Specimen Collection Kit has been specifically developed for the collection of "known" DNA from the underage mother through buccal swabs and the surgically aborted embryonic or fetal tissue. This kit allows law enforcement officials to identify the biological father through DNA.

Item.	Description
FS-DNAT/B0GEN	Criminal Paternity DNA Specimen Collection Kit, 10 kits/case

FS-DNATB0GEN	
KIT COMPONENTS	
1 ea.	Pre-sealed Kit Box (8 ⁵ / ₈ " L x 5" W x 2 ⁵ / ₈ " H) with chain of custody
1 ea.	For Law Enforcement Official Envelope affixed to bottom of kit box containing the Law Enforcement Official's Instruction Sheet (printed on yellow stock)
1 ea.	Medical Personnel Instruction Sheet/FDA Insert (printed front and back on blue stock)
1 ea.	Criminal Paternity Specimen Collection Kit Information Form printed on 3-ply NCR paper
1 ea.	Victim Buccal Swab Information Card with a swab envelope affixed
1 ea.	Pair of Nitrile Barrier Gloves
1 ea.	Sterile, cotton-tipped Swabs (2/pkg)
1 ea.	Police Evidence Seal for sealing swab envelope
1 ea.	Disposable Fingerprint Inking Sheet
1 ea.	Ziplock Bag containing a moisture-absorbing desiccant packet with biohazard label affixed
1 ea.	100 ml Tissue Specimen Bottle with embryonic/fetal tissue specimen label
1 ea.	Leakproof, tamper-evident, self-sealing Specimen Bag containing a liquid-absorbing sheet
1 ea.	Tamper-evident Kit Box Shipping Seal
1 ea.	Ziplock Bag used to protect kit box from condensation when kit is shipped to the laboratory

Insulated Specimen Shipping Box for FS-DNATB0GEN



Our 13³/₄" L x 10¹/₂" W x 13¹/₂" H reusable Insulated Specimen Shipping Box (inside dimensions: 10¹/₂" L x 7³/₄" W x 10" H) is constructed of 275 lb. test fiberboard with 1¹/₂" thick styrofoam interior. When the box is used in conjunction with refrigerant gel packs or dry ice, specimens will remain refrigerated or frozen for several days.

The Specimen Shipping Box is provided with two refrigerant gel pack labels affixed to the outside of the box and two "cover-up" dry ice labels which are affixed over the gel pack labels when dry ice is used to keep the specimen frozen while in transit to the laboratory.

Item.	Description
FS-COOLERSTY(T/B)	Insulated Specimen Shipping Box

Tip

Our boxes comply with overnight carrier regulations for shipping refrigerated or frozen biological specimens.

Refrigerant Gel Packs

These 8" L x 5¹/₂" W x 1¹/₄" H Gel Packs are safe, nontoxic and reusable. When used with our Insulated Specimen Shipping Box described above, they will help keep the collected tissue refrigerated for several days, therefore preventing degradation of the tissue specimens. We recommend six to seven gel packs be used when shipping the Criminal Paternity DNA Specimen Collection Kit to the laboratory

Item.	Description
FS-ICEPACK24	Refrigerant Gel Packs, 12 packs/case



DNA Bloodstain Specimen Collection Kit (and/or other body fluids)

Several years ago a high-profile case, which drew the attention of the media, was lost because of improper bloodstain collection procedures and poor chain of custody of the samples collected at the crime scenes.

Our Bloodstain Collection Kit is designed to eliminate these problems. Each kit is supplied with a detailed, step-by-step instruction sheet which eliminates confusion at the time of collection, ensures proper chain of custody, and provides the crime laboratory with properly collected, preserved, and clearly identifiable samples for analysis.



Tip For simplified shipping of the samples to the crime laboratory, a box version (6" L x 4¾" W x 2½" H) of our DNA Bloodstain Specimen Collection Kit is also available as Reorder No.: BSC-0GEN(B)

FS-BSC0GEN & FS-BSC0GEN(B)

KIT COMPONENTS			
1 ea.	Tamper-evident, pre-sealed Kit Envelope, 6" x 9"		
1 ea.	Kit Instruction Sheet/FDA Insert		
1 ea.	Pair of Nitrile Barrier Gloves		
1 ea.	Distilled Water Dispenser for collection of dry blood samples		
1 ea.	Police Evidence Seal for resealing kit envelope after specimen collection		
5 ea.	Sterile, cotton-tipped 3" Swabs (2/pkg)		
5 ea.	Swab Boxes		
5 ea.	Swab Box Envelopes to prevent cross-contamination		
		Item.	Description
		FS-BSC0GEN	DNA Bloodstain Specimen Collection Kit, 25 kits/case
		BSC0GEN(B)	DNA Bloodstain Specimen Collection (Box) Kit, 25 kits/case

MISSING PERSON'S FAMILY DNA COLLECTION KIT

Our Missing Person's Family DNA Collection Kit is used in cases where a family member has been reported missing.

With our kit the investigating officer collects reference DNA samples from the missing person's biological parent(s), sibling, or children, and then submits the kit to the crime laboratory. The lab then enters the results in the family member database in the Missing Person's Index of the National DNA Database (CODIS). Having the biological family member(s) tested in advance expedites the identification of recovered remains.



FS-DNA0MPF

KIT COMPONENTS	
1 ea.	Ziplock Bag, 5" x 8", with 2" x 1½" biohazard label affixed and containing a desiccant packet
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Missing Person's Family DNA Collection Kit Folder, 8½" x 6¾", with swab envelope affixed to page 3 of folder
1 ea.	Pair of Nitrile Barrier Gloves
1 ea.	Package of sterile, cotton-tipped Swabs (2/pkg)
1 ea.	Glassine Paper, 6" x 6", folded to fit in swab envelope
1 ea.	Tamper-evident Security Seal for sealing swab envelope

Tip This kit can also be used for DNA identification of remains in mass casualty situations, such as plane crashes or train wrecks.

Item.	Description
FS-DNA0MPF	Missing Person's Family DNA Collection Kit, 25 kits/case

CHILD IDENTIFICATION KIT

Child safety is a concern of parents, communities, and police throughout the country. When a child is reported missing, law enforcement agencies require accurate and current information about the child as quickly as possible in order to begin a search. In most cases, however, parents are so distraught they have difficulty answering even simple questions concerning their child.

Our Child Identification Kit is designed to eliminate these problems by providing police with photographs, thumbprints, and other pertinent information to assist in a rapid recovery of a missing or abducted child. Our kit also collects DNA saliva specimens which, in a worst-case scenario, would be used to make a positive identification of the child's remains.

Our kit can be used in conjunction with police or civic public relations programs at schools, county fairs, etc. These kits can be custom printed with your organization's name and logo on the back of each kit envelope, along with your corporate sponsor's names and logos.

Our Child Identification Kit is retained by the parent/guardian for safekeeping and only given to the police if the child is reported missing or abducted.



Tip

Quantity discounts are available. Please call for more information.



FS-CI0GEN(S)

KIT COMPONENTS


1 ea.	Kit Envelope, 6" x 9"
1 ea.	Child Identification Form, 11" x 17"
1 ea.	Pre-sealed Saliva Swabs Envelope (permanently affixed to page 4 of form)
1 ea.	Tamper-evident Security Seal for resealing saliva envelope
1 ea.	Biohazard Label for saliva swabs envelope
1 ea.	Sterile, cotton-tipped Swabs (2/pkg)
1 ea.	Pre-inked Fingerprint Sheet
1 ea.	FDA Product Insert

Item.	Description
FS-CI0GEN(S)	Child Identification Kit, 25 kits/case

FINGERNAIL SCRAPING, SWABBING, AND CLIPPINGS EVIDENCE COLLECTION KITS

Over the past few years, there has been a renewed interest in the collection of foreign material found under the fingernails of both victims and assailants. This resurgence has been generated by the ever growing use of DNA evidence by more and more investigating officers throughout the world.

The kits described below have been specifically developed to assure the proper collection, preservation, and documentation of this valuable evidence.

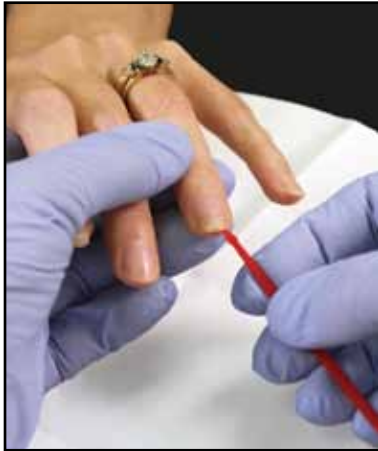

Scraping Procedure



Swabbing Procedure




Clipping Procedure

Fingernail Scraping Evidence Collection Kit

Our Fingernail Scraping Evidence Collection Kit is used in cases where the victim scratched the assailant's skin during the assault and visible debris is noted under the victim's fingernails. After collection the sealed kit is sent to the DNA section of your crime laboratory for comparison with the assailant's "known" DNA.

FS-FINSCOGEN	
KIT COMPONENTS	
1 ea.	Tamper-evident, pre-sealed Kit Envelope, 5½" x 7½"
1 ea.	Kit Instruction Sheet
1 ea.	Right Hand Envelope containing one paper bindle and one plastic fingernail scraper
1 ea.	Left Hand Envelope containing one paper bindle and one plastic fingernail scraper
1 ea.	Tamper-evident Police Evidence Seal for resealing kit envelope after specimen collection

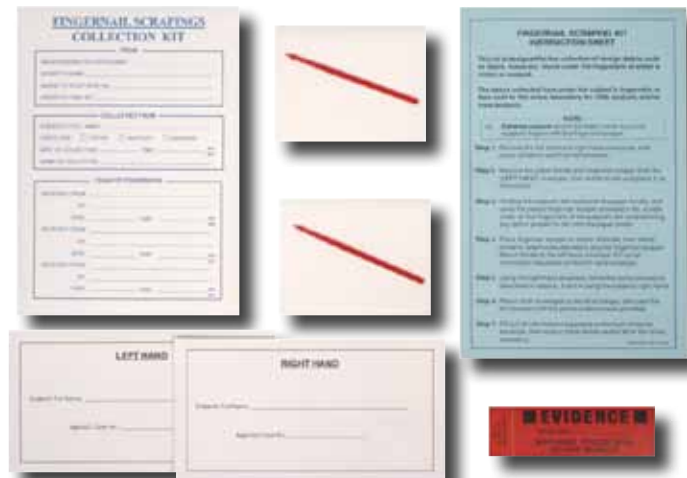

Item.	Description
FS-FINSCOGEN	Fingernail Scraping Evidence Collection Kit, 25 kits/case

Fingernail and Fingertip Swabbing Evidence Collection Kit

Our Fingernail and Fingertip Swabbing Evidence Collection Kit is used in cases when the victim states he/she scratched the assailant's skin but there is no visible debris noted under the victim's fingernails.

Simply moisten the swabs with the distilled water provided in the kit, then thoroughly swab under the victim's nails. The moistened swab will collect any trace amounts of the assailant's DNA. When the kit is submitted to the crime laboratory, the collected DNA is extracted from the swab bulbs and compared to the assailant's "known" DNA.

FS-FINSWÖGEN	
KIT COMPONENTS	
1 ea.	Tamper-evident, pre-sealed Kit Envelope, 5½" x 7½"
1 ea.	Kit Instruction Sheet
2 ea.	Sterile, cotton-tipped Swabs (2/pkg)
1 ea.	3 ml Distilled Water
1 ea.	Right Hand Swab Box
1 ea.	Left Hand Swab Box
1 ea.	Tamper-evident Police Evidence Seal for resealing kit envelope after specimen collection




Item.	Description
FS-FINSWÖGEN	Fingernail and Fingertip Swabbing Evidence Collection Kit, 25 kits/case

Fingernail Clippings Evidence Collection Kit

During an assault the victim's fingernail(s) may break off during the struggle. Our Fingernail Clippings Evidence Collection Kit is used to collect the remaining fingernail for comparison with the broken fingernail found where the assault took place.

When the clipping(s) is submitted to the crime laboratory, trace amounts of the assailant's DNA may be found on the clipping and on the recovered broken fingernail as well.

FS-FINCLOGEN	
KIT COMPONENTS	
1 ea.	Tamper-evident, pre-sealed Kit Envelope, 5½" x 7½"
1 ea.	Kit Instruction Sheet
1 ea.	Metal Fingernail Clipper
1 ea.	Right Hand Envelope containing one paper bindle
1 ea.	Left Hand Envelope containing one paper bindle
1 ea.	Tamper-evident Police Evidence Seal for resealing kit envelope after specimen collection

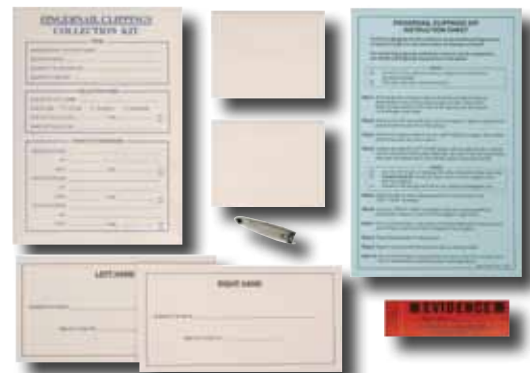

Item.	Description
FS-FINC-0GEN	Fingernail Clippings Evidence Collection Kit, 25 kits/case

DNA FTA® MATRIX CARDS

Many of our custom kit customers are now using FTA® cards in their **I.D. TECHNOLOGIES** DNA CODIS Database Kits and Sexual Assault Evidence Collection Kits for the collection and preservation of whole blood and buccal cells/saliva specimens.

I.D. TECHNOLOGIES is now providing the FTA® cards and the accessories, described below, directly to crime laboratories and law enforcement agencies. FTA® technology is a revolutionary tool for the simplified collection, purification, archiving, and analysis of pure DNA and RNA in whole blood and buccal cells/saliva.


FTA® cards are composed of a sophisticated filtration matrix impregnated with a patented formulation of powerful protein denaturants, a chelating agent, and a free radical trap designed to protect and entrap nucleic acids. It is nontoxic to humans yet prevents the growth of bacteria and other microorganisms. FTA® Matrix Cards also rapidly inactivate organisms such as blood borne pathogens, and, therefore, offer user safety. FTA® technology maintains the integrity of the nucleic acids contained in the specimen from nucleases, oxidation, and UV damage.

FTA® Matrix Cards permit the collected specimens to be stored at room temperature without degradation for up to 14 years of documented data (as of 2005). This unique feature enables the user to purify nucleic acids from samples at their convenience. In addition to being an excellent archival storage media for nucleic acids, products containing FTA® technology can be transported easily and tracked conveniently through the U.S. Postal Service without the need for temperature control or hazardous/restricted packaging.

FTA® Elute Cards

•	Provides stable, long-term room temperature DNA storage with easy release using simple elution method
•	Eliminates lengthy and multiple step isolation procedure
•	Indicator card changes from purple to white when the card comes in contact with the sample.
•	FTA® Elute Card Overall dimensions: 1.5 x 3 inches with four 11 mm. diameter circles on a 1.5 x 1.5 inch matrix; each circle holds 40 µL.
•	FTA® Indicating Elute Card Overall dimensions: 1.5 x 3 inches with a single one-inch diameter circle on a 1.5 x 1.5 inch matrix; circle holds 40 µL.


Item.	Description
FS-BODPAFTAELU	FTA® Elute Micro Cards, 100/pkg
FS-BODPAFTAELU-I	FTA® Indicating Elute Micro Cards, 100/pkg

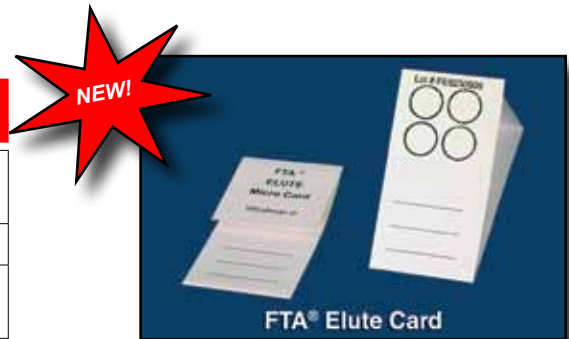

EasiCollect™ Devices

•	Unique design engineered for reliable, uniform collection and application of cells to the included FTA® Card
•	Foam collection surface captures large number of cells for analysis
•	EasiCollect™ device holds FTA® card securely in place while molded pressure clips ensure constant pressure for reproducible transfer of cells
•	Compatible with manual or automated punching systems

Item.	Description
FS-BODPAFTAEEASI	EasiCollect™ Devices, 50 cards/pkg

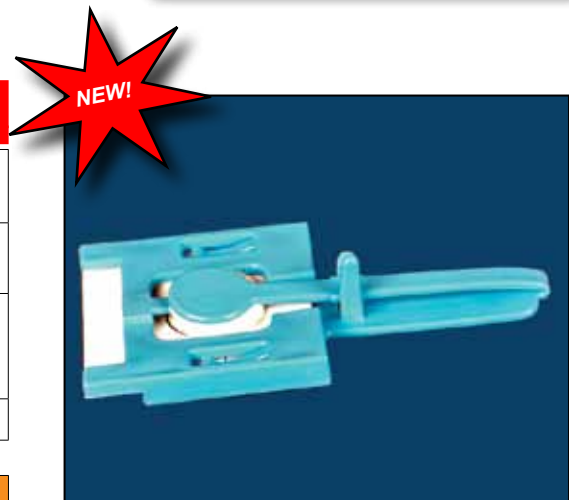


FTA® is a registered trademark of Whatman Inc.

FTA® Cards for whole blood and buccal cells/saliva

•	Designed for the collection, transportation, storage, and purification of nucleic acids from whole blood and buccal cells/saliva
•	Integral cover sheet protects the sample and permits the card to stand on its own while drying
•	FTA® <u>Classic</u> Card Overall dimensions: 3 x 5 inches with four one-inch diameter circles on a 3 x 3 inch matrix; each circle holds 100-125 µL
	FTA® <u>Mini</u> Card Overall dimensions: 1.88 x 2.88 inches with two one-inch diameter circles on a 1.25 x 2.88 inch matrix; each circle holds 100-125 µL
	FTA® <u>Micro</u> Card Overall dimensions: 1.5 x 3 inches with a single one-inch diameter circle on a 1.5 x 1.5 inch matrix; circle holds 100-125 µL



Item.	Description
FS-BODPAFTA4	FTA® Classic Cards, 100 cards/pkg
FS-BODPAFTA2	FTA® Mini Cards, 100 cards/pkg
FS-BODPAFTA1	FTA® Micro Cards, 100 cards/pkg

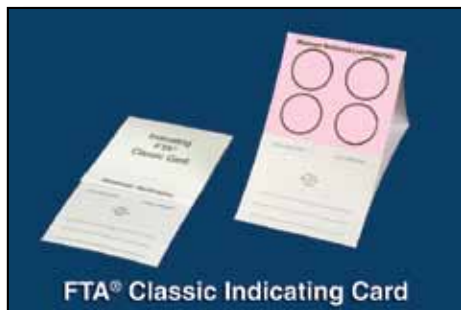

FTA® Indicating Cards for buccal cells/saliva

•	Indicator card changes from pink to white when the card comes in contact with the sample. The dye washes out during processing and does not inhibit analysis or transformation
•	Integral cover sheet protects the sample and permits the card to stand on its own while drying
•	FTA® <u>Classic Indicating</u> Card Overall dimensions: 3 x 5 inches with four one-inch diameter circles on a 3 x 3 inch matrix; each circle holds 100-125 µL
	FTA® <u>Mini Indicating</u> Card Overall dimensions: 1.88 x 2.88 inches with two one-inch diameter circles on a 1.25 x 2.88 inch matrix; each circle holds 100-125 µL
	FTA® <u>Micro Indicating</u> Card Overall dimensions: 1.5 x 3 inches with a single one-inch diameter circle on a 1.5 x 1.5 inch matrix; circle holds 100-125 µL

Item.	Description
FS-BODPAFTA4IND	FTA® Classic Indicating Cards, 100 cards/pkg
FS-BODPAFTA2IND	FTA® Mini Indicating Cards, 100 cards/pkg
FS-BODPAFTA1IND	FTA® Micro Indicating Cards, 100 cards/pkg

FTA® is a registered trademark of Whatman Inc.

ADDITIONAL COMPONENTS

Harris Uni-Core™ Punch and Mats

- Designed for the quick and easy sampling of dried materials contained within FTA® Card
- Razor sharp cutting tip has a hollow body with a plunger for expelling the punched sample
- Sharp cutting edge ensures clean, precise, uniform sampling
- Precision-tooled tip: polished steel, rust-resistant, heat-treated, and case hardened to Rockwell 60; when handled properly, tip will provide 1500 punches when used with a cutting mat
- Can be sterilized with ethylene oxide or autoclaved for 20 minutes at 250° F/15 psi
- Tip can be rinsed in ethanol after each use to eliminate cross-contamination
- 4.125 inches in overall length
- Available in 1.25 mm, 2.0 mm, 3.0 mm, and 6.0 mm diameters
- Cutting mats (two sizes available) are 3 mm thick. Self-healing inert mat ensures crisp sample cuts and protects and extends the life of the cutting edge of the tool (can be sterilized with ethanol)




Item.	Description
FS-PUNCH1.25	1.25 mm Uni-Core™ Punch, 4 punches/pkg
FS-PUNCH2	2.0 mm Uni-Core™ Punch, 4 punches/pkg
FS-PUNCH3	3.0 mm Uni-Core™ Punch, 4 punches/pkg
FS-PUNCH6	6.0 mm Uni-Core™ Punch, 4 punches/pkg
FS-MAT2.5X3	Small Cutting Mat, 2½" x 3", 1 punch/pkg
FS-MAT6X8	Large Cutting Mat, 6" x 8", 1 punch/pkg

Fingerstick Lancets



- Haemolance® 21 gauge, fixed depth, sterile lancets are used in all our kits that call for a fingerstick blood sample. They are quick, simple to use, and safe
- To obtain a blood sample, simply twist off the protective base of the lancet; then hold the Haemolance® against the subject's finger and press the green button
- The 21 gauge stick immediately retracts into the body of the lancet when the button is released guarding against accidental contact with the stick
- After use discard the Haemolance® into a Sharps container

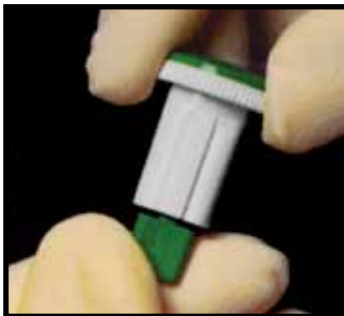

Item.	Description
FS-LANCET21H	21 Gauge Lancet, 150 lancets/box

Other sizes are also available. Call for more information

903 Protein Saver™ Card

- The sample collection area of the 903 Protein Saver™ Card contains five ½" circles
- Each circle holds 75 to 80 µl of sample
- Overall card dimensions: 3½" X 2⅛"




Item.	Description
FS-BODPAPS	903 Protein Saver™ Card, 100 cards/pkg

BUCCAL CELLS/SALIVA COLLECTION SWABS

Sterile Omni-Swabs

- Noninvasive device for collection of saliva and buccal cells
- Unique brushlike swab head composed of a nonreactive, absorbent filter paper easily ejects from the stem of the swab for transfer of samples into tubes and multi-well plates
- Sterilized and individually wrapped

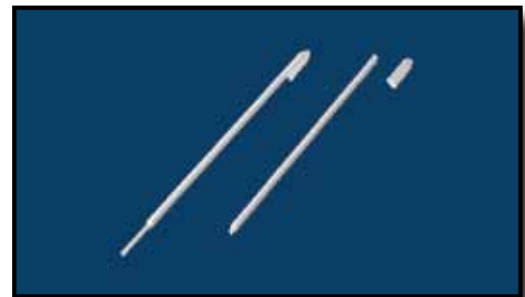

Item.	Description
FS-SWABOMNI	Omni-Swabs, 100 swabs/box

Sterile Foam-Tipped Swabs


- A nonabrasive device for the collection of saliva and buccal cells to be transferred to FTA® Cards
- Plastic swab shaft measures 5" long; foam tip measures 2⁹/₃₂" wide by 1¹/₈" length
- Sterilized and individually wrapped



Item.	Description
FS-SWABFOAM	Foam-Tipped Swabs, 50 swabs/box

Cap-Shure™ Sterile Swabs

- Sterilized and individually wrapped
- Cotton-tipped with either a plastic or wood shaft
- Reclosable, ventilated cap to protect the sample while allowing it to air-dry

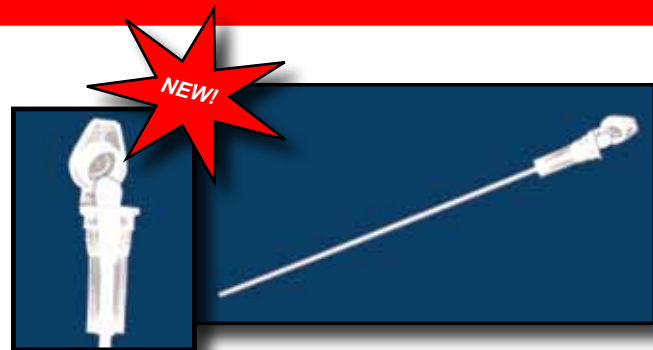

Item.	Description
FS-SWABCAPP	Cap-Shure™ Sterile Swabs with Plastic Shaft, 50 per pack (shown)
FS-SWABCAPW	Cap-Shure™ Sterile Swabs with Wood Shaft, 50 per pack

FTA® CARD MAILING AND STORAGE

FTA® Card Mailer

- For use with the Multi-barrier Pouches to transport FTA® Matrix Cards in accordance with U.S. Postal Service
- Dimensions: 4.5 x 7.125 inches (11.4 x 18.1 cm)


Item.	Description
FS-MAILERFTA	FTA® Card Mailer, 50 mailers/pkg



FTA® Shipping & Storage Multi-barrier Pouches

- For storing FTA® Classic and Classic Indicating Cards
- Dimensions:
 - Large (for all Classic FTA® cards) - 4.375 x 6.5 inches (11.11 x 16.51 cm)
 - Small (for all smaller FTA® cards) - 3.75 x 3 inches (9.53 x 7.62 cm)
- Constructed of seven laminated layers that protect the card from exposure to any gaseous or liquid external element
- Tamper-evident seal maintains sample security
- Outer paper layer for labeling or writing

Item.	Description
FS-POUCH4375X65	Large Multi-barrier Pouch, 100 pouches/pkg
FS-POUCH375X3	Small Multi-barrier Pouch, 100 pouches/pkg

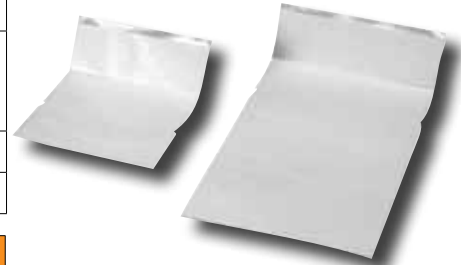

Indicating Desiccant Packets

- Ensures the FTA® matrix remains dry during storage
- Measures .75" W x 1.75" L x .25" H
- One gram cobaltous silica gel with millions of pores to soak up water vapor and eliminate moisture and dampness
- Changes to pink upon absorption of moisture
- Complies with FDA guidelines

Item.	Description
FS-DESICCANTIND	Indicating Desiccant Packets, 1000 desiccants/case



Desiccant Packets

- Measures .5" W x 1" L x .25" H
- .75 gram silica gel with millions of pores to soak up water vapor and eliminate moisture and dampness

Item.	Description
FS-DESICCANT.75	Desiccant Packets, 1000 desiccants/case



BLOOD AND URINE SPECIMEN COLLECTION KITS (FOR ALCOHOL AND/OR DRUG DETERMINATIONS)



The series of DWI Specimen Collection Kits designed and manufactured by I.D. TECHNOLOGIES are supplied with a detailed step-by-step instruction sheet which eliminates confusion at the time of specimen collection by clearly defining and separating the responsibilities of the investigating officer from those of the medical personnel collecting the specimens.

Our kits assure crime laboratory personnel that they will receive undamaged, properly collected, and clearly identifiable specimens for analysis. If you mail your kits to the crime laboratory, our boxed kits comply with commercial carrier and U.S. Postal regulations pertaining to the shipment of biological specimens.

We realize that some agencies have unique requirements for their DWI kits. I.D. TECHNOLOGIES has the expertise to design kits that meet your exact requirements, whether you need pre-addressed kit boxes or the inclusion of your own evidence submission forms. Give us a call – we'll help you select the appropriate components and provide you with a detailed price and delivery quotation for your review.

Urine Specimen Collection Kits

FS-UC1 (One Bottle)

KIT COMPONENTS

1 ea.	Pre-sealed Kit Box (5" L x 3¾" W x 2½" H) with biohazard label affixed
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Specimen Bottle Holder (polyurethane foam)
1 ea.	100 ml Urine Specimen Bottle
1 ea.	Ziplock Bag containing a liquid-absorbing sheet
1 ea.	Investigating Officer's Report and Chain of Possession Label for the ziplock bag
1 ea.	Tamper-evident Specimen Security/Identification Seal
1 ea.	Tamper-evident Kit Box Shipping Seal for resealing kit box after specimen collection



FS-UC2 (Two Bottles)

KIT COMPONENTS

1 ea.	Pre-sealed Kit Box (6" L x 4¾" W x 2½" H) with biohazard label affixed
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Specimen Bottle Holder (polyurethane foam)
2 ea.	100 ml Urine Specimen Bottle
1 ea.	Ziplock Bag containing a liquid-absorbing sheet
1 ea.	Investigating Officer's Report and Chain of Possession Label for the ziplock bag
2 ea.	Tamper-evident Specimen Security/Identification Seals
1 ea.	Tamper-evident Kit Box Shipping Seal for resealing kit box after specimen collection



Item.	Description
FS-UC1	Urine Specimen Collection Kit (One Bottle), 25 kits/case
FS-UC2	Urine Specimen Collection Kit (Two Bottles), 25 kits/case

Tip

Our FS-UC1 and FS-UC2 Urine Specimen Collection Kits are specifically designed for DWI specimen collection; however, we also provide similar kits for pre-employment and routine drug testing of law enforcement and governmental personnel. In those kits a urine temperature monitoring thermometer is affixed to the bottle, eliminating the need for direct observation at the time of collection. Please contact our Customer Service Representatives at 800.438.7884 for more information about our Workplace Substance Abuse Products.

Blood Specimen Collection Kit

FS-BA2 (Two Blood Tubes)

KIT COMPONENTS

1 ea.	Pre-sealed Kit Box (5" L x 3¾" W x 2½" H) with biohazard label affixed
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Subject Consent Form
1 ea.	Inner Box Blood Tube Holder with Investigating Officer's Report and Chain of Possession printed on security flap
1 ea.	Tamper-evident Evidence Seal for sealing inner box
2 ea.	10 ml Blood Collection Tubes (gray-stoppered) containing 100 mg of sodium fluoride and 20 mg of potassium oxalate
1 ea.	20 g x 1½" Multi-sample Needle
1 ea.	Safety Tube and Needle Holder
1 ea.	PVP Prep Pad (alcohol-free)
1 ea.	Ziplock Bag containing a liquid-absorbing sheet
2 ea.	Tamper-evident Specimen Security/Identification Seals
1 ea.	Tamper-evident Kit Box Shipping Seal for resealing outer kit box after specimen collection



Blood and Urine Specimen Collection Kit

FS-BA3 (2 Blood Tubes & 1 Bottle)

KIT COMPONENTS

1 ea.	Pre-sealed Kit Box (6" L x 4¾" W x 2½" H) with biohazard label affixed
1 ea.	Kit Instruction Sheet/FDA Insert
1 ea.	Subject Consent Form
1 ea.	Inner Box Blood Tube and Urine Bottle Holder with Investigating Officer's Report and Chain of Possession printed on security flap
1 ea.	Tamper-evident Evidence Seal for sealing inner box
1 ea.	100 ml Urine Specimen Bottle
2 ea.	10 ml Blood Collection Tubes (gray-stoppered) containing 100 mg of sodium fluoride and 20 mg of potassium oxalate
1 ea.	20 g x 1½" Multi-sample Needle
1 ea.	Safety Tube and Needle Holder
1 ea.	PVP Prep Pad (alcohol-free)
1 ea.	Ziplock Bag containing a liquid-absorbing sheet
3 ea.	Tamper-evident Specimen Security/Identification Seals
1 ea.	Tamper-evident Kit Box Shipping Seal for resealing outer kit box after specimen collection



Item.	Description
FS-BA2	Blood Specimen Collection Kit (Two Blood Tubes), 25 kits/case
FS-BA3	Blood and Urine Specimen Collection Kit (Two Blood Tubes and One Bottle), 25 kits/case

Urine Specimen Collection Kit

If you collect urine specimens at your facility and hand deliver them to the laboratory for testing, our UC-0DOC kit contains all the necessary components to meet your requirements.

FS-UC0DOC	
KIT COMPONENTS	
1 ea.	Product Identification Sheet/FDA Insert
1 ea.	Pre-sealed 120 ml widemouth Collection Cup/Bottle with a urine temperature monitoring thermometer affixed
1 ea.	Leakproof, tamper-evident, self-sealing, double-pouch Specimen/Document Bag, 6" x 10", containing a liquid absorbing sheet
1 ea.	Urine specimen bottle Security Seal

Item.	Description
FS-UC0DOC	Urine Specimen Collection Kit, 25 kits/case



(Each UC-0DOC is packaged in a tamper-evident, heat-sealed integrity bag.)

Self-Sealing, Leakproof Specimen Bags

Our Specimen Bags are designed to contain liquid if a urine bottle should leak. The bags utilize a unique adhesive strip which forms a leakproof, tamper-evident seal that is impossible to reopen without destroying the bag. These bags are available with and without liquid absorbing sheets which contain an immobilizing agent that "gels" when it comes in contact with urine. Unlike liquid absorbing packets, our sheets will absorb from a few drops to the entire contents of a specimen bottle.

Our Double-Pouch Specimen Bags have two compartments to separate the specimen bottle(s) from the documentation/request for analysis form.



Item.	Description
FS-SPEBAGS	Urine Specimen Collection Kit, 25 kits/case
FS-SPEBAGS/A	4" x 6" bag w/liquid absorbing sheet, 200 bags/pkg
FS-SPEBAGDRUG	6" x 10" double-pouch bag, 200 bags/pkg
FS-SPEBAGDRUG/A	6" x 10" double-pouch bag w/liquid absorbing sheet, 200 bags/pkg

Bubble Pack Blood Tube Protectors

Our self-sealing Bubble Pack Blood Tube Protectors are available in one tube and two tube versions. The two tube version has two compartments to prevent the glass tubes from touching each other. These blood tube protectors are an inexpensive way to ensure blood tubes will not break in shipping or storage.

Item.	Description
FS-BUBPA3X55	Single Blood Tube Protector, 3" X 5½", 250 protectors/case
FS-BUBPA4X55CS	Double Blood Tube Protector, 4" X 5½", 250 protectors/case



Liquid Absorbing Sheets

Our Liquid Absorbing Sheets contain an immobilizing agent that gels when it comes in contact with blood, urine, or other body fluids.

Unlike liquid absorbing packets, our 4¼" long by 2¼" wide sheets will absorb from a few drops up to 200 ml of liquid. These are the same absorbing sheets used in all of our Specimen Collection Kits which collect blood and urine.

Item.	Description
FS-ABSHEET2.25X425	Liquid Absorbing Sheets, 100 sheets/pkg



BLOOD SPECIMEN COLLECTION TUBES

The blood tubes described below are routinely used by medical personnel to collect victim and suspect blood specimens for a variety of forensic testing. All blood tubes shipped from **I.D. TECHNOLOGIES** have an expiration date of 12 months or greater.

Gray-Stoppered

Our 10 ml Gray-Stoppered Glass Blood Tubes (16 x 100 mm) contain 100 mg of sodium fluoride and 20 mg potassium oxalate. Our Gray-Stoppered Tubes are used for collecting blood specimens for alcohol and/or drug testing.



Item.	Description
FS-TUG10	Gray-Stoppered Blood Collection Tubes, 100 tubes/pkg

Lavender-Stoppered

These 7 ml Lavender-Stoppered Glass Blood Tubes (13 x 100 mm) contain a 15% solution of EDTA (K3) which is the preferred preservative for DNA typing of collected blood specimens.



Item.	Description
FS-TUL7L	Lavender-Stoppered Blood Collection Tubes, 100 tubes/pkg

Yellow-Stoppered

These 7 ml Yellow-Stoppered Glass Blood Tubes (13 x 100 mm) contain ACD Solution B* which is the preferred preservative for blood typing.



Item.	Description
FS-TUY7	Yellow-Stoppered Blood Collection Tubes, 100 tubes/pkg

*ACD Solution B of trisodium citrate, 13.2 g/L; citric acid, 4.8 g/L; and dextrose 14.7 g/L, 1.0 mL

Red-Stoppered

Our 10 ml Red-Stoppered Glass Blood Tubes (16 x 100 mm) are silicone coated and contain no additives.



Item.	Description
FS-TUR10	Red-Stoppered Blood Collection Tubes, 100 tubes/pkg

Blood Tube Specimen Security Identification Seals

Our extremely Tamper-evident Blood Tube Security Seals leave a void print on the tube and seal if there is an attempt to remove the 7" long by 1" wide seal. These seals can be easily written on with ballpoint pens.



Item.	Description
FS-SV165A	Blood Tube Specimen Security Seal, 100 seals/pkg

URINE SPECIMEN COLLECTION SUPPLIES

105 ml Urine Specimen Collection Cups



Our 105 ml Urine Specimen Collection Cups are ideal for urine collection from both male and female subjects. Each cup is supplied in a heat-sealed integrity pouch and is available with or without urine temperature monitoring thermometers.

- Screenprinted with graduated volume markings
- Made of polystyrene
- Diameter of cup opening is 2½"
- Can also be used for on-site testing

Item.	Description
FS-105CT	105 ml Urine Specimen Collection Cups w/thermometer, 125 cups/ case
FS-105C	105 ml Urine Specimen Collection Cups w/o thermometer, 125 cups/case

240 ml Urine Specimen Collection Cups



These 240 ml Urine Specimen Cups come with a leak-resistant snap and lock lid, and each cup is supplied in a heat-sealed integrity pouch. They are available with or without a urine temperature monitoring thermometer affixed to the outside of the cup.

- Ideal for on-site testing
- Wide, stable base (2½")
- Widemouth design (3" dia.) facilitates urine specimen collection from female subjects

Item.	Description
FS-240CT	240 ml Urine Specimen Cups w/thermometer, 140 bottles/case
FS-240C	240 ml Urine Specimen Cups w/o thermometer, 140 bottles/case

Urine Specimen Collection and Transport Jars



Our pre-sealed widemouth (2" dia) 120 ml jar is perfect for female subjects because it eliminates the need to transfer the urine from a collection cup to a specimen transport bottle.

Item.	Description
FS-120JT	120 ml jar w/thermometer, 150 jars/case
FS-120J	120 ml jar w/o thermometer, 150 jars/case

Urine Specimen Collection and Transport Bottles



Our 60 ml and 100 ml Urine Specimen Bottles are pre-sealed at I.D. TECHNOLOGIES via an internal bottle cap integrity seal and are supplied with volume markings. Both the 60 ml and 100 ml urine specimen bottles are available with or without urine temperature monitoring thermometers affixed to the outside of the bottle.

Item.	Description
FS-60BT/M	60 ml bottle w/thermometer, 250 bottles/case
FS-60B/M	60 ml bottle w/o thermometer, 250 bottles/case
FS-100BT/M	100 ml bottle w/thermometer, 250 bottles/case
FS-100B/M	100 ml bottle w/o thermometer, 250 bottles/case

Temp-Secur® (Patent# 07/968,385)

This unique patented design prevents subjects from bringing in someone else's drug free urine and using an external heat source to bring the specimen to an acceptable temperature.

These 100 ml Urine Specimen Bottles have a temperature monitoring thermometer that is concealed from the subject's view by an external label with a tear-strip. After receiving the specimen from the subject, the collector simply pulls off the tear-strip to expose the specimen temperature. Each bottle is supplied in a heat-sealed integrity pouch.




Item.	Description
FS-100BT/TS	Temp-Secur® 100 ml Urine Specimen Bottles, 100 bottles/case

100BT/TS	Temp-Secur® 100 ml Urine Specimen Bottles, 100 bottles/case	\$116.00
----------	---	----------

Specimen Security/Identification Seals

Our urine bottle seals are tamper-evident and provide a clearly identifiable sample for the lab. Each seal is 7" x 1¼" with a 5/8" tab.

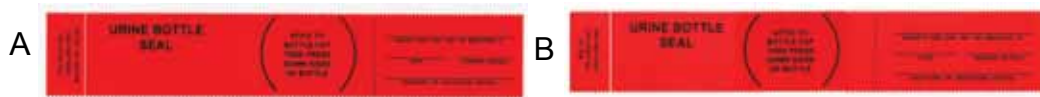

Item.	Description
FS-ST278	(A) 60 ml Specimen Security/Identification Seals, 250 seals/pkg
ST219	(B) 100 ml Specimen Security/Identification Seals, 250 seals/pkg

OnTrak TesTstik™ On-Site Urine Drug Screening Kits

If your policy calls for on-site screening, we offer a wide variety of single and multiple OnTrak TesTstik™ test sticks that combine unparalleled reliability and ease of use in one advanced, cost-effective design. Their superior latex technology and patented manufacturing process virtually eliminate false positives and provide results as accurate as laboratory gas chromatography/mass spectrometry (GC/MS) tests. There's no mixing, clock-watching, or specialized operator training needed. Best of all, they are affordable and are available in one, two, and three drug configurations for today's most commonly abused drugs. The OnTrak TesTstik™ will provide you with accurate, easy-to-read results within five minutes, and your findings will stay accurate and readable for up to 30 minutes after testing is completed.

On-Site Urine Drug Screening Kits - TEST STICKS

FEATURES

- Easy Storage - Stores at room temperature for as long as three years from date of manufacture
- Simple to Use - No mixing of reagents or pipetting required
- Accurate - Meets current legally defensible SAMHSA* cutoff detection levels
- Clear - Easy-to-interpret test results appear next to the name of the drug you are testing for (blue band indicates a negative result)
- Long Life - Three year expiration date from date manufactured

* Substance Abuse Mental Health Service Administration



Item.	Description
FS-OTTSAMP	OnTrak TesTstik™ for Amphetamines, 50 kits/case
FS-OTTSBAR	OnTrak TesTstik™ for Barbiturates, 50 kits/case
FS-OTTSBENZ	OnTrak TesTstik™ for Benzodiazepines, 50 kits/case
FS-OTTSCOC	OnTrak TesTstik™ for Cocaine, 50 kits/case
OTSMETH	OnTrak TesTstik™ for Methamphetamines, 50 kits/case
OTTSOPI	OnTrak TesTstik™ for Morphine 300 ng/mL (Opiates), 50 kits/case
OTTSPCP	OnTrak TesTstik™ for PCP, 50 kits/case
OTTSTHC	OnTrak TesTstik™ for THC 50 ng/mL (Marijuana), 50 kits/case
OTTS2	OnTrak TesTstik™ 2 for Cocaine and THC, 50 kits/case
OTTS3	OnTrak TesTstik™ 3 for Cocaine, THC, and Morphine 300 ng/mL (Opiates), 50 kits/case

OnTrak TesTstik™ is a registered trademark of Varian, Inc.